


Experiment : 5

Experiment :- To determine the amount of Fe^{2+} and Fe^{3+} ions by permanganometry.

Apparatus - Pipette, burette, beakers, conical flask, burette stand and clamp.

Chemicals - Mohr's salt solution (Ferrous ammonium sulphate, $\text{FeSO}_4 (\text{NH}_4)_2 \text{SO}_4 \cdot 6\text{H}_2\text{O}$), permanganate (KMnO_4) and sulphuric acid (H_2SO_4)

Chemical Equation -

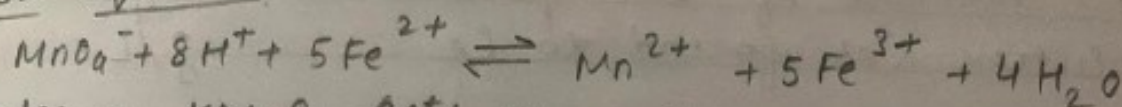

Indicator - KMnO_4 acts as self indicator.

Observation

a) Standardization of KMnO_4 solution.

Volume of 0.1 N FAS (N_1) solution taken for each titration = 10 ml.

S.No.	Burette reading (ml)		Vol. of KMnO_4 used 1 ml
	Initial	Final	
1	0	10.5	10.5
2	0	10.5	10.5
3	0	10.5	10.5

Mean Volume (V_1) = 10.5 ml

b) Determination of Fe^{2+}

Volume of given sample taken for each titration = 10 ml

S.No.	Burette reading (ml)		Vol. of KMnO_4 used (ml)
	Initial	Final	
1	0	10.5	10.5
2	0	10.5	10.5
3	0	10.5	10.5

Khushi Gopal

Experiment: To determine the amount of Fe^{2+} and Fe^{3+} ions by permanganometry.

Apparatus: pipette, burette, beakers, conical flask, burette stand and clamp.

Chemicals: Mohr's Salt solution (Ferrous Ammonium sulphate; $\text{FeSO}_4 (\text{NH}_4)_2 \text{SO}_4 \cdot 6\text{H}_2\text{O}$), permanganate (KMnO_4) and sulphuric acid (H_2SO_4)

Theory: Mn^{7+} oxidises Fe^{2+} in acidic medium to Fe^{3+} and itself gets reduced to divalent chromium (Mn^{2+})
$$\text{MnO}_4^- + 8\text{H}^+ + 5\text{Fe}^{2+} \rightleftharpoons \text{Mn}^{2+} + 5\text{Fe}^{3+} + 4\text{H}_2\text{O}$$

KMnO_4 acts as a self indicator. If Fe^{3+} is present in the original solution, it can be reduced by boiling the soln with zinc pieces in acidic medium and can be titrated with standard KMnO_4 . The endpoint in this case corresponds to presence of both Fe^{2+} and Fe^{3+} ions in the solution.

PROCEDURE. i) Standardisation of KMnO_4

- 1) Transfer 10ml of the standard 0.1N ferrous ammonium sulphate (FAS) solution to a clean conical flask using a pipette.
- 2) Add 5ml of 4N sulphuric acid.
- 3) Titrate the solution against KMnO_4 solution taken in a burette. The color of the solution changes from colorless to pink.
- 4) Note the volume of the solution used and repeat the titration atleast 5 times and take the mean of closely related readings (4).

Teacher's Signature: Kushigod

Mean volume (V_2) = _____ ml

(c) Determination of Fe^{2+} and Fe^{3+} (Total iron content)
 volume of given sample taken for each titration
 = 10 ml

S.No	Burette reading (ml)		Volume of $KMnO_4$ used (ml)
	Initial	Final	
1.	0	14.5	14.5
2.	0	14.5	14.5
3.	0	14.5	14.5

Mean Volume (V_3) = 14.5 ml

CALCULATIONS

i) Normality of $KMnO_4$ solution

$$(KMnO_4) \quad N_1 V_1 = N_2 V_2 \text{ (FAS)} = 0.1 N \times 10 \text{ ml}$$

$$\text{Normality of } KMnO_4 (N_1) = \frac{0.1 \times 10}{V_1} = S(N) = \frac{1}{10.5} N = 0.095 N$$

ii) Determination of Fe^{2+}

Volume of solution taken = 10 ml

Volume of $KMnO_4$ solution used (V_2) = 10.5 ml

Normality of $KMnO_4$ (N_1) = 0.095 N

$$\text{Normality of } Fe^{2+} = N = \frac{V_2 \times S}{V_1} = \frac{10.5 \times 0.095}{10} \approx 0.1$$

$$\text{Strength of } Fe^{2+} = \frac{56}{10} \times N = 56 \times \frac{0.1}{10} = 5.6 \text{ gm/L}$$

iii) Determination of Fe^{+3} in a mix of Fe^{2+} and Fe^{3+}

Volume of sol. taken used = 10 ml

Volume of $KMnO_4$ used (V_3) = 14.5 ml

Normality of $KMnO_4$ (N_1) = 0.095 N

$$\text{Normality of total } Fe^{+3} \quad N_1 = \frac{V_3 \times S}{10} = \frac{0.095 \times 14.5}{10} \approx 0.14 N$$

$$\text{Strength of total } Fe = N_1 \times \text{Eq wt.} = 56 \times 0.14 = 7.7 \text{ gm/L}$$

Khushi goel

i) Determination of Fe^{2+} .

- 1) Pipette out 10 ml of the given substance in a titration conical flask.
- 2) Add 5 ml of 4 N sulphuric acid.
- 3) Titrate the soln. against KMnO_4 soln taken in burette. The color of the substance changes from colorless to pink.
- 4) Note the volume of the soln used and repeat the titration atleast 5 times and take the mean of closely related reading (V_2).

ii) Determination of Fe^{2+} and Fe^{3+} (Total iron content)

- 1) Pipette out 10 ml of aq. soln into the conical flask [The given solution has already been boiled with 2-3 g of zinc pieces and 5 ml of dilute H_2SO_4 to reduce Fe^{3+} to Fe^{2+}].
- 2) Add 5 ml of 4 N H_2SO_4 .
- 3) Titrate it with standard KMnO_4 solution till the color changes from colorless to pink.
- 4) Note the volume of the solution used and repeat the titration atleast 5 times and take the mean of the closely related readings (V_3).

General Calculations

i) Normality of KMnO_4 soln

Applying the normality solution

$$N_1 V_1 (\text{KMnO}_4) = N_2 V (\text{FAS}) = 0.1 \text{ N} \times 10 \text{ ml.}$$

$$\text{Normality of } \text{KMnO}_4 \quad N_1 = \frac{0.1 \times 10}{V_1} = 5 (\text{N})$$

ii) Determination of Fe^{2+}

Vol of solution taken = 10 ml.

Teacher's Signature :

Khushi
Gautam

$$= N_1 \times 56 = \text{gm/L} = 7.7 \text{ gm/L}$$

Strength of Fe^{3+} ions = Normality of Eq wt

$$\frac{(V_3 - V_2) \times 8 \times 56}{10} = \underline{2.1} \text{ gm/L}$$

Results : The amount of $\text{Fe}^{2+} = \underline{5.6} \text{ gm/L}$
 The amount of $\text{Fe}^{3+} = \underline{2.1} \text{ gm/L}$

Khushi
 400

Vol. of KMnO_4 solution used = V_2 ml

Normality of KMnO_4 = 8 N

Normality of Fe^{2+} $N_1 = \frac{V_2 \times 8}{10}$

Strength = $N_1 \times \text{Eq. wt} \Rightarrow N_1 \times 56 \text{ gm/L}$

iii) Determination of Fe^{3+} in a mixture of Fe^{2+} and Fe^{3+}

Volume of soln taken = 10 ml

Volume of KMnO_4 solution used = V_3 ml

Normality of KMnO_4 = 8 N

Normality of total Fe = $N_1 = V_3 \times 8 / 10$

Strength of total Fe = $N_1 \times \text{Eq. wt} = N_1 \times 56 \text{ gm/L}$

Strength of Fe^{3+} ions = Normality \times Eq. wt
 $(V_3 - V_2) \times 8 \times 56 / 10$

Results \rightarrow The amount of Fe^{2+} = 5.6 gm/L
 The amount of Fe^{3+} = 2.1 gm/L

Expected CLO's / Daily life Application: In quantitative^{to} and qualitative determination of Fe^{2+} and/or Fe^{3+} present in an ore or compound, water sample, etc.