


Experiment - Determine  $pK_a$  value of acetic acid by pH metric titration.

Apparatus - Pipette, burette, beakers, funnel, pH meter and glass electrode.

Chemicals Required - Sodium hydroxide ( $NaOH$ ) and acetic acid ( $CH_3COOH$ )

Chemical Equations -  $HA_{(aq)} + NaOH_{(aq)} \rightarrow Na^+ A_{(aq)} + H_2O$

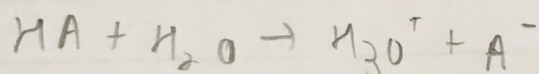

Observations - Normality of standard  $NaOH = 0.1N$


Vol. of $NaOH$ added from the burette (mL)		pH
1	0.5	2.85
2	0.7	3.77
3	1.2	3.96
4	1.7	4.10
5	2.2	4.24
6	2.7	4.36
7	3.2	4.44
8	3.7	4.56
9	4.2	4.65
10	4.7	4.75
11	5.2	4.85
12	5.7	4.92




Experiment-8

Aim - Determine  $pK_a$  value of acetic acid by pH meter titration.

Theory - A pH meter will be used to follow the titration of an unknown weak acid,  $HA(aq)$  with sodium hydroxide,  $NaOH(aq)$

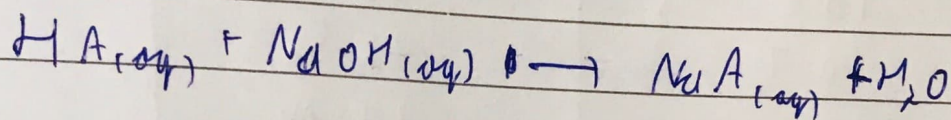

The weak acid has a concentration around 0.1 M. The result of the pH versus volume of  $NaOH$  plot is "S" shaped curve which is not as steep as the one arising from the titration of strong acid. The equivalence point (the time) will be at alkaline pH (not 7 as in strong vs strong base). From the equivalence point, the concentration of an unknown acid  $HA$  is found. In addition, the acid constant  $K_a$  can be determined.

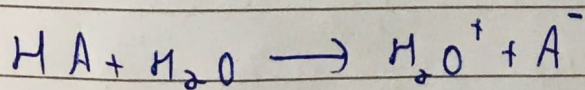

$$\left\{ \begin{array}{l} pH = pK_a + \log \frac{[salt form]}{[acid]} \\ \text{Henderson-Hasselbalch eq} \end{array} \right.$$

$$K_a = \frac{[H_3O^+][A^-]}{[HA]}$$

Procedure - A). Titration of unknown  $HA$  with standard  $NaOH$ :

1. Calibrate the pH meter with standard buffer sol<sup>n</sup> of pH=4

Teacher's Signature \_\_\_\_\_



13.	6.2	5.05
14.	6.7	5.18
15.	7.2	5.31
16.	7.7	5.46
17.	8.2	5.62
18.	8.7	5.97
19.	9.2	9.86
20.	9.7	12.01
21.	10.2	12.36
22.	10.7	12.52
23.	11.2	12.64
24.	11.7	12.70
25.	12.2	12.76
26.	12.7	12.82
27.	13.2	12.87
28.	13.7	12.90
29.	14.2	12.94
30.	14.7	12.95

Draw a graph pH vs volume of NaOH (sample graph, as shown). Find out pKa value of  $\text{CH}_3\text{COOH}$  from the graph as under:-

NaOH (ml) at equivalence point	NaOH (ml) at half equivalence point	pH at half equivalence point
$x$ (along x-axis) 9.1	$\frac{x}{2}$ (along x-axis) $\frac{9.1}{2} = 4.55$	pH at $\frac{x}{2}$ = A (along y-axis) 4.8



of 9, then rinse the glass electrode and immerse it in the beaker. Position the beaker so that the titrant can be easily added.


2. Pipette out 50 ml of acid and into a clean beaker, dip the glass electrode. Record the pH.

3. Initially, add 0.5 ml of 0.1 NaOH solution at a time, record the pH change is more than 0.2-0.3 units, then start adding 0.2 ml of NaOH each time (i.e. near to the equivalence point, decrease the volume of NaOH added) so that the change in pH is small enough to yield a good shape of plot.

4. After the rapid change in pH (after the equivalence point), the volume of NaOH may again be increased to 0.5 ml per addition. Make at least 10 more additions after the equivalence point so that the region with the plateau can be plotted.

5.  $pK_a$  is determined by examining the titration curve. The negative log of  $K_a$  is  $pK_a$  and is same as the pH at half the volume of equivalence point. ( $pH = pK_a$  when logarithm term is zero which is then is zero and  $[salt] = [acid]$ ). This is true at half, equivalence point i.e. Henderson Hasselbalch equation.



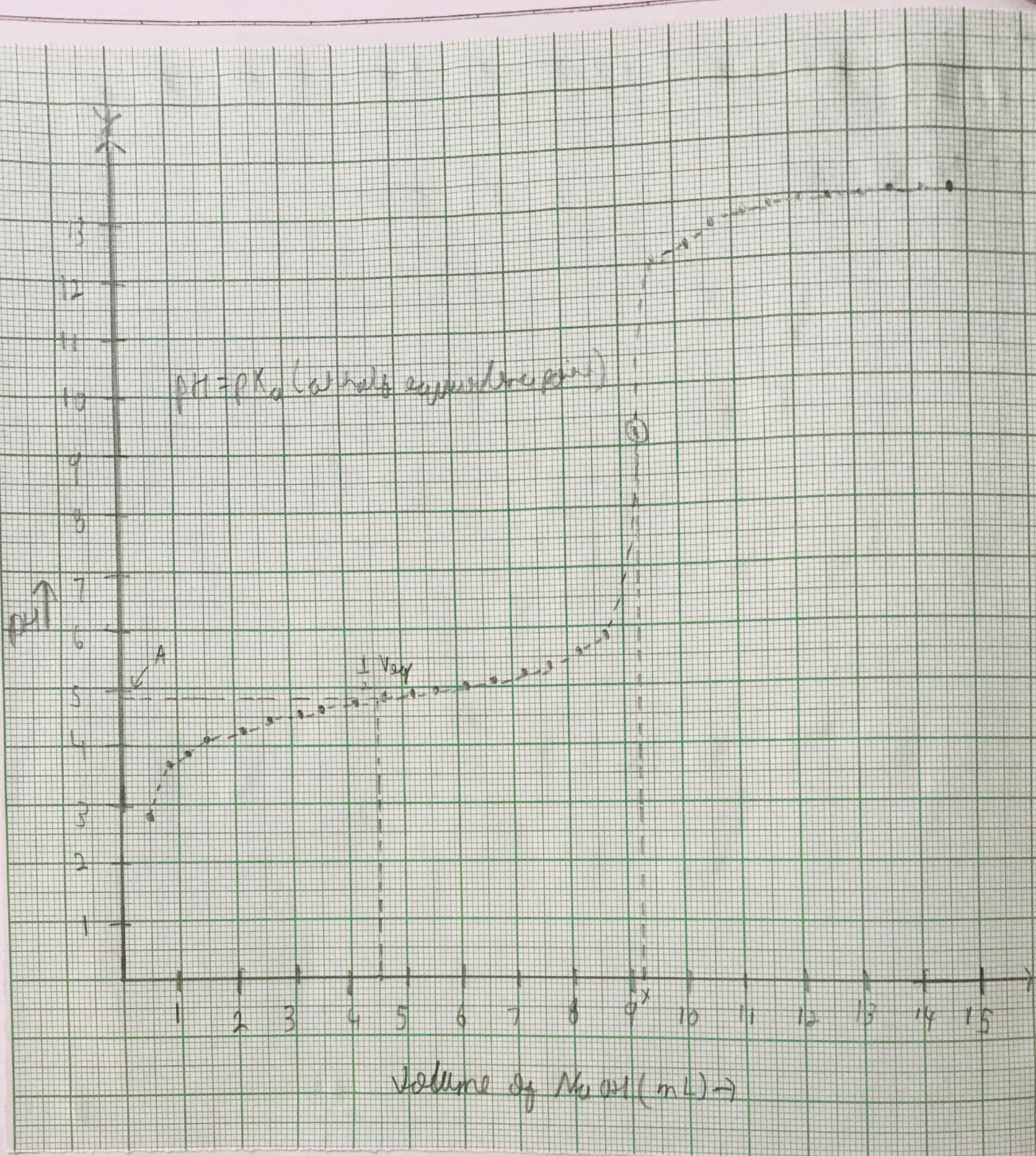

Result- The  $\text{pK}_a$  of acetic acid is 4.8



Result - The  $pK_a$  of acetic acid is 4.8.

Precautions - (i) Rinse the pipette/burette with the sol<sup>n</sup> to be transferred to the titration flask/burette

(ii) Do not rinse the titration flask.