

Date : 23<sup>rd</sup> nov, 2020


Experiment No:- 08

Experiment :- Determine  $pK_a$  value of acetic acid by pH metric titration.

Apparatus :- Pipette, burette, beakers, funnel, burette stand, clamp, pH meter and glass electrode.

Chemicals Required :- Sodium hydroxide (NaOH) and acetic acid ( $CH_3COOH$ )

Chemical Equations:-

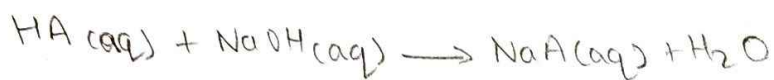


Experiment - 8→ Aim:-

Determine  $pK_a$  value of acetic acid by pH metric titration.

→ Theory:-

A pH meter will be used to follow the titration of an unknown weak acid,  $HA(aq)$  with sodium hydroxide,  $NaOH(aq)$

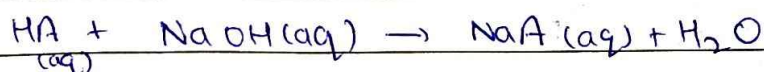

The weak acid has a concentration around 0.1 M. The result of the pH versus volume of NaOH plot is "S" shaped curve which is not as steep as the one arising from the titration of strong acid. The equivalence point (this time) will be at alkaline pH (not 7 as in strong acid vs strong base). From the equivalence point, the concentration of an ~~unkn~~ unknown acid HA is found. In addition, the acid constant  $K_a$  can be determined.



$$pH = pK_a + \log \frac{[\text{salt form}]}{[\text{acid form}]}$$

Henderson-Hasselbalch equation

$$K_a = \frac{[H_3O^+][A^-]}{[HA]}$$

Teacher's Signature \_\_\_\_\_

(Nishay)

Observations:-

Normality of standard  $\text{NaOH} = 0.1\text{N}$

Volume of $\text{NaOH}$ added from the burette (mL)		pH
1.	0.4	2.85
2.	0.9	3.77
3.	1.4	3.96
4.	1.9	4.10
5.	2.4	4.24
6.	2.9	4.36
7.	3.4	4.4
8.	3.9	4.56
9.	4.4	4.65
10.	4.9	4.75
11.	5.4	4.85
12.	5.9	4.92
13.	6.4	5.05
14.	6.9	5.18
15.	7.4	5.31
16.	7.9	5.46
17.	8.4	5.62
18.	8.9	5.97
19.	9.4	9.36
20.	9.9	12.01
21.	10.4	12.36
22.	10.9	12.52
23.	11.4	12.64
24.	11.9	12.70

Arthay



→ Procedure:-# Titration of ~~known~~ unknown HA with standard NaOH:-

- 1) Calibrate the pH meter with the standard buffer solution of  $\text{pH} = 4$  or  $9$ , then rinse the glass electrode and immerse it in the beaker. Position the burette so that the titrant can be easily added.
- 2) Pipette out  $50\text{ ml}$  of acetic acid into a clean beaker, dip the glass electrode. Record the pH.
- 3) Initially, add  $0.5\text{ ml}$  of  $0.1\text{ N}$  NaOH solution at a time, record the pH (after each titration) until the pH change is more than  $0.2 - 0.3$  units, then start add  $0.2\text{ ml}$  of NaOH each time (i.e., near to the equivalence point, decrease the volume of NaOH added) so that the change in pH is small enough to yield a good shape of plot.
- 4) After the rapid change in pH (after the equivalent point), the volume of NaOH may again be increased to  $0.5\text{ ml}$  per addition. Make at least 10 more additions after the equivalence point so that the region with the plateau can be plotted.
- 5)  $\text{pK}_a$  is determined by examining the titration curve. The negative log of  $K_a$  is  $\text{pK}_a$  and is same as the pH at half the volume of equivalence point.  
( $\text{pH} = \text{pK}_a$  when logarithm term is zero which in term turn is zero once  $[\text{salt}] = [\text{acid}]$ . This is true at half equivalence point) cf. Henderson - Hasselbalch

Teacher's Signature \_\_\_\_\_

Nishu

Volume of NaOH added from the burette (mL)		pH
25.	12.4	12.76
26.	12.9	12.82
27.	13.4	12.87
28.	13.9	12.90
29.	14.4	12.94
30.	14.9	12.95

Draw a graph pH vs. volumes of NaOH (sample graph is as shown). Find out  $pK_a$  value of  $CH_3COOH$  from the graph as under :-

NaOH (mL) at equivalence point	NaOH (mL) at half equivalence point	pH (at half equivalence point) = $pK_a$
9.9	4.4	4.65



equation.

→ Result:-

The  $pK_a$  of acetic acid is 4.65.


→ Precautions:-

- (i) Rinse the pipette / burette with the solution to be transferred to the titration flask / burette.
- (ii) Do not rinse the titration flask.

Teacher's Signature \_\_\_\_\_

*Arjun*



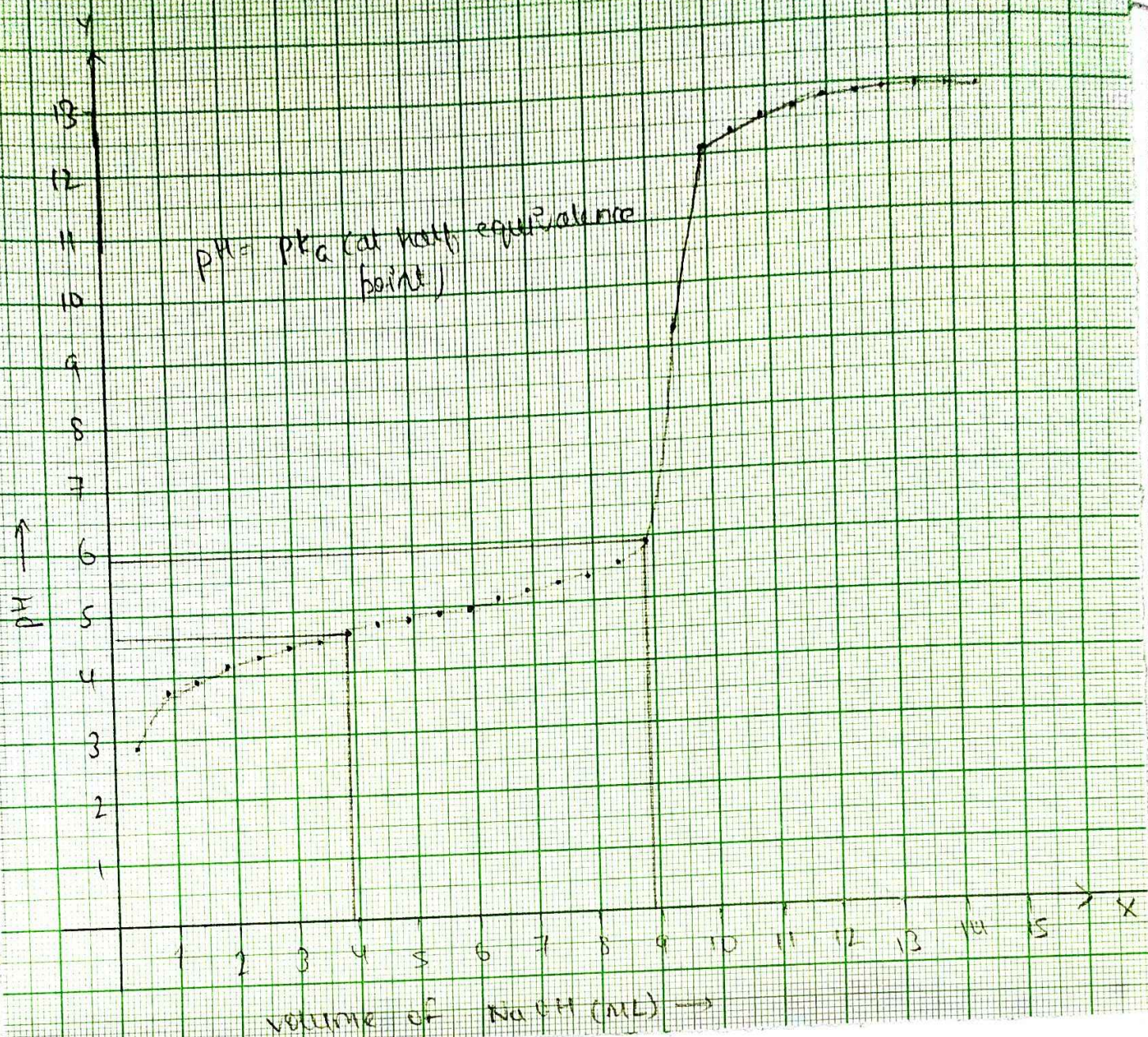

Result :-

The  $\text{pKa}$  of acetic acid is 4.65.