

Experiment - 5

Experiment - to determine the amount of Fe^{2+} and Fe^{3+} ions by permanganometry.

Apparatus - Pipette, burette, beakers, conical flask, burette stand and clamp.

Chemicals - Mohr's salt solⁿ, $\text{FeSO}_4(\text{NH}_4)_2\text{SO}_4 \cdot 6\text{H}_2\text{O}$, permanganate (KMnO_4) & Sulphuric acid (H_2SO_4).

Chemical Equations - $\text{MnO}_4^- + 8\text{H}^+ + 5\text{Fe}^{2+} \rightleftharpoons \text{Mn}^{2+} + 5\text{Fe}^{3+} + 4\text{H}_2\text{O}$

Indicator - KMnO_4 acts as self indicator.

Observation - 1. Standardization of KMnO_4 solⁿ.

Volume of 0.1N $\text{FeSO}_4(\text{NH}_4)_2\text{SO}_4$ solⁿ taken for each titration = 10 ml

S. No.	Burette reading (ml)		Vol of KMnO_4 used (ml)
	Initial	Final	
1.	0	10.3	10.3
2	0	10.3	10.3
3	0	10.3	10.3


Mean volume (V_1) = 10.5 ml

Experiment: - To determine the amount of Fe^{2+} & Fe^{3+} ions by permanganometry.

Apparatus: Pipette, burette, flasks, conical flask, burette stand & clamps

Chemicals: Mohr's salt solⁿ (Ferrous ammonium sulphate) $FeSO_4 (NH_4)_2 SO_4 \cdot 6H_2O$, Permanganate ($KMnO_4$) and sulphuric acid H_2SO_4

Theory: Mn^{7+} oxidises Fe^{2+} in acidic medium to Fe^{3+} and itself gets reduced to divalent chromium (Mn^{2+})

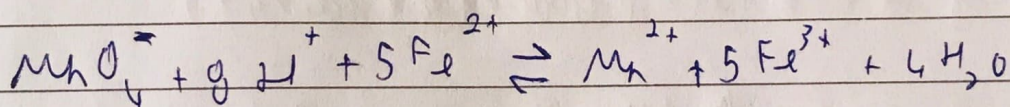

$KMnO_4$ acts as a self indicator. If Fe^{3+} is present in the original solⁿ, it can be reduced by boiling the soln with zinc pieces in acidic medium and can be titrated with standard $KMnO_4$. The endpoint in this case corresponds to presence of both Fe^{2+} and Fe^{3+} ions in the solⁿ.

Procedure: 1. Standardisation of $KMnO_4$ -

1. Transfer 10 ml of the standard 0.1 N ferrous ammonium sulphate (FAS) solⁿ to a clean conical flask using a pipette.

2. Determination of Fe^{2+}

Vol. of given sample taken for each titration = 10 ml

S.No	Burette Reading (ml)		Volume of KMnO_4 used (ml)
	Initial	Final	
1	0	10.3	10.3
2	0	10.3	10.3
3	0	10.3	10.3

Mean vol (V_2) = 10.3 ml

3. Determination of Fe^{2+} & Fe^{3+} (total iron content)

volume of given sample taken for each titration = 10 ml.

S.No	Burette Reading (ml)		Vol. of KMnO_4 used (ml)
	Initial	Final	
1	0	14.3	14.3
2	0	14.3	14.3
3	0	14.3	14.3

Mean vol (V_3) = 14.3 ml

(ii) Add 5 ml of 4 N sulphuric acid.

(iii) Titrate the solⁿ against KMnO_4 solⁿ taken in a burette the colour of the solⁿ changes from colourless to pink.

(iv) Note the volume of the solⁿ used and repeat the titration atleast 5 times and take the mean of closely related readings (V_1).

2. Determination of Fe^{2+}

(i) Pipette out 10 ml of the given substance in a titration conical flask.

(ii) Add 5 ml of 4 N sulphuric acid

(iii) ~~Titrate~~ Titrate the solⁿ against KMnO_4 solⁿ taken in burette the colour of the substance changes from colourless to pink.

(iv) Note the volume of the solⁿ used and repeat the titration atleast 5 times and take the mean of closely related reading (V_2).

3. Determination of Fe^{2+} and Fe^{3+} (Total iron content)

(i) Pipette ~~out~~ out 10 ml of aq. solⁿ into the conical

Calculations:- 1. Normality of KMnO_4 solⁿ

$$(\text{KMnO}_4) N_1 V_1 = N_2 V_2 \text{ (FAS)} = 0.1 \text{ N} \times 10 \text{ ml}$$

$$\text{Normality of } \text{KMnO}_4 (N_1) = \frac{0.1 \times 10}{V_1} = \frac{0.1 \times 10}{10.3} = S(N) = 0.097 \text{ N}$$

2. Determination of Fe^{2+}

Vol of solⁿ taken = 10 ml

Vol of KMnO_4 solⁿ used (V_2) = 10.3 ml

Normality of $\text{KMnO}_4 (N_1) = 0.097 \text{ N}$

$$\text{Normality of } \text{Fe}^{2+} = N = V_2 \times S = 10.3 \times 0.097 \approx 0.1$$

$$\text{Strength of } \text{Fe}^{2+} = 56 \times N = 56 \times 0.1 = 5.6 \text{ gm/L}$$

3. Determination of Fe^{+3} in a mixed Fe^{2+} & Fe^{3+}

Vol of solⁿ taken used = 10 ml

Vol. of KMnO_4 used (V_3) = 14.3 ml

Normality of $\text{KMnO}_4 (N_1) = 0.097 \text{ N}$

$$\text{Normality of total Fe} \mid N_1 = \frac{V_3 \times S}{10} = \frac{0.097 \times 14.3}{10} \approx 0.138$$

$$\text{Strength of total Fe} = N_1 \times \text{Eq. wt} = 0.138 \times 56 = 7.728 \text{ m/L}$$

flask (The given solⁿ has already been boiled with 2-3 g of zinc pieces and 5 ml of dilute H_2SO_4 to reduce the Fe^{3+} to Fe^{2+}).

(ii) Add 5 ml of 4N H_2SO_4 .

(iii) Titrate it with standard $KMnO_4$ solⁿ till the colour changes from colourless to pink.

(iv) Note the volume of the solⁿ used and repeat the titration at least 5 times and take the mean of the closely related readings (V_3).

General calculation:

1. Normality of $KMnO_4$ solⁿ

Applying the normality solⁿ

$$N_1 V_1 (KMnO_4) = N_2 V_2 (FAS) = 0.1 N \times 10 \text{ ml}$$

$$\text{Normality of } KMnO_4 \quad N_1 = \frac{0.1 \times 10}{V_1} = 5 (N)$$

2. Determination of Fe^{2+}

Vol of solⁿ taken = 10 ml.

Strength of Fe^{3+} ions = Normality of Eq wt

$$\frac{(V_3 - V_2) \times N \times 56}{10} = 2.1 \text{ gm/L}$$

Result-

The amount of $Fe^{2+} = 5.6 \text{ gm/L}$

The amount of $Fe^{3+} = 2.1 \text{ gm/L}$

Vol. of $KMnO_4$ soln used = V_2 ml

Normality of $KMnO_4$ = $8N$

$$\text{Normality of } Fe^{2+} \quad N_1 = \frac{V_2 \times 8}{10}$$

$$\text{Strength} = N_1 \times \text{Eq. wt} = N_1 \times 56 \text{ gm/L}$$

3. Determination of Fe^{3+} in a mixture of Fe^{2+} and Fe^{3+}

Volume of soln taken = 10 ml

Volume of $KMnO_4$ soln used = V_3 ml

Normality of $KMnO_4$ = $8N$

$$\text{Normality of total Fe} = N_1 = \frac{V_3 \times 8}{10}$$

$$\text{Strength of total Fe} = N_1 \times \text{Eq. wt} = N_1 \times 56 \text{ gm/L}$$

$$\text{Strength of } Fe^{3+} \text{ ius} = \text{Normality} \times \text{Eq. wt.}$$

$$= (V_3 - V_2) \times 8 \times \frac{56}{10}$$

Result - The amount of Fe^{2+} = 5.6 gm/L

The amount of Fe^{3+} = 2.1 gm/L